




2 A Teschio e denti, 1489.

RL 19058 v (B. 41 r; K/P 42 v).

Il vacuo della cassa dell'occhio, e 'r vacuo dell'osso sostenitore della guancia, e quello del naso e della bocca sono d'eguale profondità e terminano sotto il senso comune [sede centrale delle sensazioni] per linea perpendicolare. E ciascuna d'esse vacuità ha tanto di profondità quant'è la terza parte del volto dell'omo, cioè dal mento a capegli.

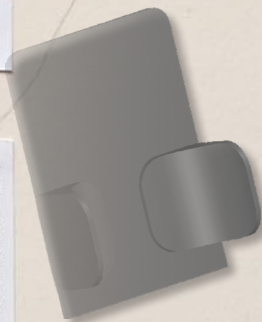

A NEW SOLUTION

TTS TRUNK TOTAL SUPPORT

Positioning kit for the trunk with modular inserts



Why Vip ?

- WITHOUT MOUNTING HARDWARE AND ALUMINIUM SHELL
- WITH TENSIONING STRAP SYSTEM
- LIGHTWEIGHT, COMFORTABLE, SAFE
- MAXIMUM SUPPORT, CONTROL AND CONTAINMENT OF THE TRUNK



STRUCTURAL KIT WITH
A FLAT BASE OF CONSTRUCTION
AND 2 INSERTS




POLYURETHANE
MEMORY FOAM PADDING



MICROCLIMATIC COVER

CLINICAL INFORMATION/DESTINATION OF USE

- Non-ambulatory users with moderate/high trunk instability that require support
- Users with limited trunk control

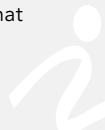

TTS EMBRACING
SHAPES

TECHNICAL FEATURES

- Structural kit composed by a flat base of construction and 2 polyethylene closed celled inserts, modifiable and configurable, useful to individually support or correct the posture of the trunk during the life time use
- Polyurethane memory foam padding
- Microclimatic ventilated cover, fireproof and antibacterial
- Tensioning strap system with H 10 cm straps for wheelchair installation

SPARE
PARTS

- Flat base of construction
- Padding

- Cover
- Tensionable system

WARNINGS

The trunk Positioning System **Inserto Back Vip**, in its design integrity and numerical and dimensional completeness of the supplied components (Structural Kit, Padding, Cover) can be easily adapted to sizes/morphology/deformities of the user. This kind of operations make the commissioning referable to a **serial device**.


Alternatively

The **commissioning** of the Positioning System **Inserto Back Vip**, deprived or modified, even partially, of its design integrity of numerical and dimensional completeness of the supplied components, built as per written prescription of a professional operator in function of the anatomy/morphology/deformities of the user, through the measurement detection and direct trials, can be referable to a **custom device**.

For partially sighted users this document can be visualised in PDF format on the web site www.versainserto.com with the possibility to enlarge the text



VIP

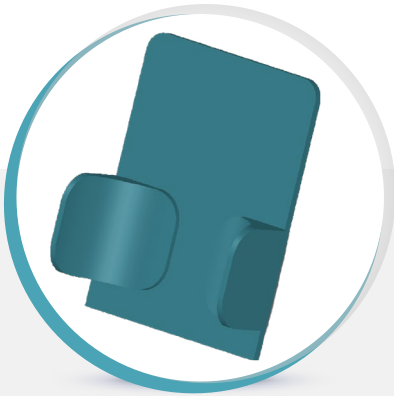

Inserto Vip

2 INSERTS MODIFIABLE AND CONFIGURABLE




MICROCLIMATIC VENTILATED COVER

° Each postural unit is ordered based upon the size of mobility/wheelchair device.

| INSERTO VIP SIZES TABLE | | | | | | | | | | SIZES (cm) | |
|--|-------|-------|-------|-------|-------|-------|-------|-------|-------|------------|-------|
| MODEL | 3235 | 3438 | 3642 | 3845 | 4045 | 4245 | 4250 | 4545 | 4550 | 4850 | 5050 |
| Effective shell width | 27 | 29 | 31 | 33 | 35 | 37 | 37 | 40 | 40 | 43 | 45 |
| Effective shell height | 35 | 38 | 42 | 45 | 45 | 45 | 50 | 45 | 50 | 50 | 50 |
| Effective wheelchair external tubes width (cm) min/max | 31/35 | 33/37 | 35/39 | 37/41 | 39/43 | 41/45 | 41/45 | 44/48 | 44/48 | 47/51 | 49/53 |

Semi-rigid lumbar containment/contouring insert

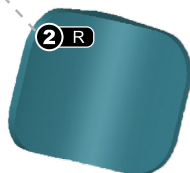


Semi-rigid lumbar containment/contouring insert



1

Rigid flat base of construction

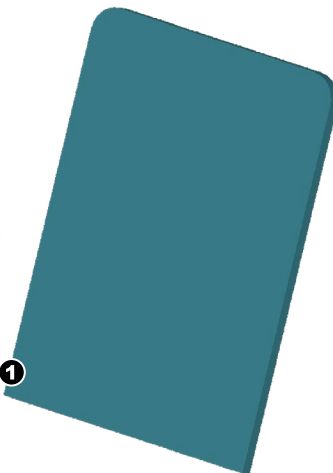



MANUFACTURED BY



Via Montagna Z.I. - Lotto 41
72023 Mesagne (BR) - ITALY
Tel. +(39) 0831 777840
Fax. +(39) 0831 730739
Email: sales@promedicare.it
www.versainserto.com
www.promedicare.eu

CERTIFIED COMPANY WITH QUALITY
MANAGEMENT SYSTEM ISO 9001 / EN ISO 13485



ISO 9001
ISO 13485
BUREAU VERITAS
Certification



DISTRIBUTED BY



CE 2 YEAR WARRANTY

www.versainserto.com