

Thinking

BASICO
inserto[®]
SEAT


SEAT



2 A Teschio e denti, 1489.

R.I. 19058 v (B. 41 r; K/P 42 v).

Il vacuo della cassa dell'occhio, e 'i vacuo dell'osso sostenitore della guancia, e quello del naso e della bocca sono d'eguale profondità e terminano sotto il senso comune [sede centrale delle sensazioni] per linea perpendicolare. E ciascuna d'esse vacuità ha tanto di profondità quant'è la terza parte del volto dell'omo, cioè dal mento a capegli.

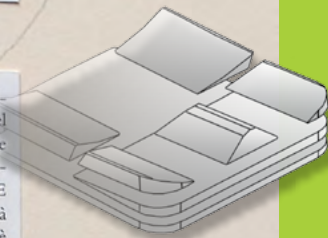

A NEW SOLUTION
INDIVIDUALLY CUSTOMISABLE

A NEW SOLUTION INDIVIDUALLY CUSTOMISABLE

Construction kit of the individually customisable seat Perfect reproduction of anatomical shapes




Why Basico ?

- LIGHT POSITIONING
- WIDER SUPPORT
- HIGH COMFORT



STRUCTURAL KIT
MODULAR AND CONFIGURABLE

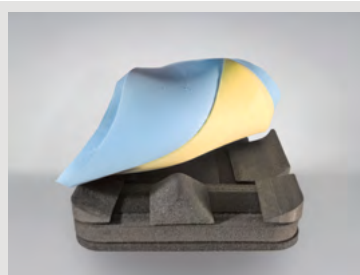

PADDING SYSTEM
Free Shape



EXCELLENT SUPPORT
AND STABILITY

CLINICAL INFORMATION / DESTINATION OF USE

The pelvis positioning system Inserto Seat Basico is recommended preferably, but not exhaustively, for users with mobility limits , in particular:

- Those who need a postural support as well as frequent positioning changes;
- Those who need a good support and a lateral, anterior containmnet of the pelvis with a medium/high risk of decubitus issues;
- Those who have leg lenght discrepancies, pelvic rotation and obliquity, posterior and anterior pelvis.

The pelvis positioning system has to be used with any supporting/wheelchair base properly fixed through suitable gripping systems.

TECHNICAL FEATURES

Each positioning system Inserto BASICO includes:

- A structural kit composed by a flat base of construction and 6 closed cell polyethylene inserts, modular and configurable, useful to support or correct the posture of the pelvis during the life time use;
- A padding, washable with a gentle detergent, which is composed by different layers of foam:
 - posterior part: two layers of foam (memory foam and flexible open celled foam with low spring back action) in order to provide the maximum comfort and the highest protection from the risk of decubitus sores;
 - anterior part: a layer of open celled foam with low spring back action to enhance the natural position of the thighs;
- New 3D cover, air-exchange, incontinent, not flammable, antibacterial. It can be washed in tepid water with a gentle detergent. Do not tumble dry. The fabric doesn't require the use of disinfectants.

SPARE PARTS

- Flat base of construction
- Cover
- Kit of padding

WARNINGS

The Positioning System Inserto Seat Basico, in its design integrity and numerical and dimensional completeness of the supplied components (structural kit, padding kit, cover) can be easily adapted to sizes/morphology/deformities of the user. This kind of operations make the commissioning referable to a serial device.

Alternatively

The commissioning of the Positioning System Inserto Seat Basico, deprived or modified, even partially, of its design integrity of numerical and dimensional completeness of the supplied components, built as per written prescription of a professional operator in function of the anatomy/ morphology/deformities of the user, through the measurement detection and direct trials, can be referable to a custom device.

For partially sighted users this document can be visualised in PDF format on the web site www.versainser.to.com with the possibility to enlarge the text.





Inserto Basico

6 MODULAR AND CONFIGURABLE INSERTS

"NEW 3D" COVER
STRECH, INCONTINENT, AIR-EXCHANGE.

° Each postural unit is ordered based upon the size of mobility/wheelchair device


INSERTO BASICO SIZES TABLE	PAEDIATRIC						ADULT						SIZES (cm)	
	UUS	USS	US	XXS	XS	XS1S	XS1	S	M	M1S	M1	L	XL	
Effective width (cm)	25	30	30	36	36	40	40	42	42	45	45	48	48	
Max width achieved with adaptation (cm)	27	32	32	38	38	42	42	44	44	47	47	50	50	
Min. width achieved with removal of material (cm)	23	28	28	34	34	38	38	40	40	43	43	46	46	
Effective length (cm)	30	30	38	38	42	40	45	45	50	45	50	45	50	
Max. length achieved with adaptation (cm)	32,5	32,5	40,5	40,5	44,5	42,5	47,5	47,5	52,5	47,5	52,5	47,5	52,5	
Min. length achieved with rear notches (cm)*	24	24	32	32	36	34	39	39	44	39	44	39	44	
Min. and max. range of anterior height of flat base of construction+ leveling insert 4 (cm)**	3,5 / 4,5		** height referred to models UUS and XL											
Posterior height of flat base of construction 5 (cm)**	1,5		** all sizes											
Weight of positioning system min/max: 0,8 kg/2,1 kg														

* possible and further reductions due to removal of material

° Max load (referred to model XL 48cm x 50cm): 135 kg

CONVERSION TABLE

INCHES	10	12	12,5	13	13,5	14	15	15,5	16	16	16,5	16,5	17	17,5	17,5	17,5	18	18,5	18,5	19	19,5	19,5	20	20	22
CM.	25	30	32	33	34	36	38	39	40	40,5	42	42,5	43	44	44,5	45	46	47	47,5	48	50	50,5	52	52,5	55

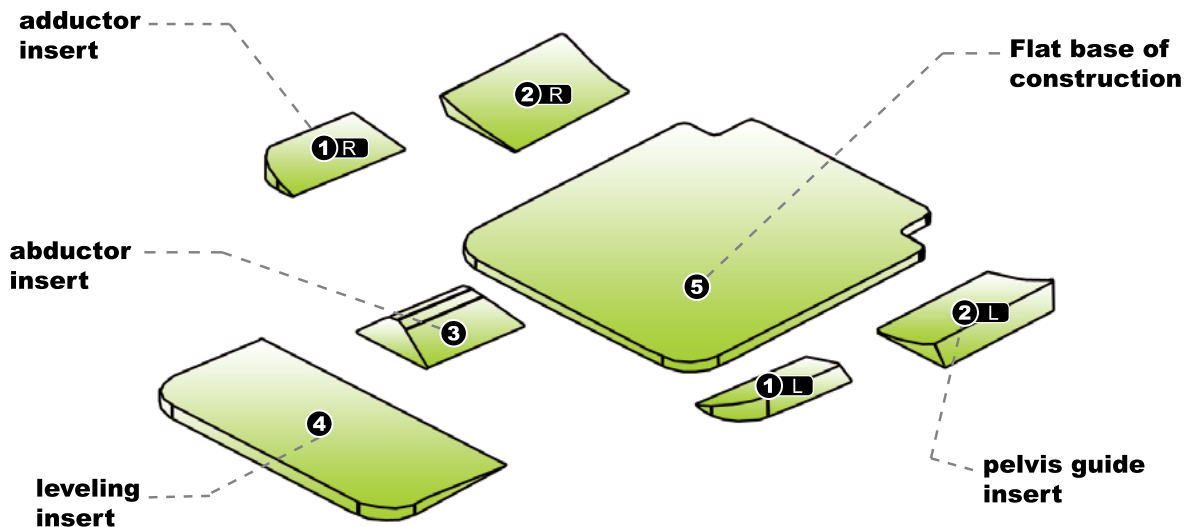


Customisations of Basico



ADJUSTMENT OF TROCHANTERIC
SUPPORT PREVENTING ISCHIAL LOADS



FIXED
ADDUCTION

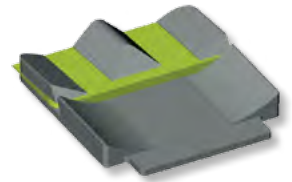


ANTI-THRUST
OF THE PELVIS




FIXED
ABDUCTION




WIND-SWEPT

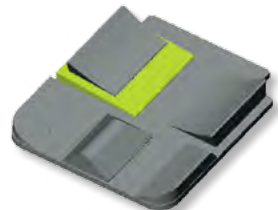

PELVIC
OBLIQUITY

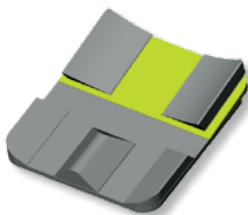

FIXED ANGLE

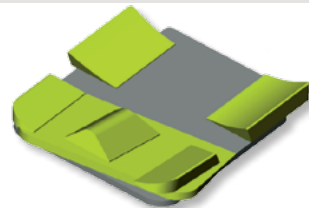

PELVIC ROTATION

INDIVIDUALLY CUSTOMISABLE

Construction kit of the individually customisable seat



THE INNOVATIVE SOLUTION

A range developed to answer the individual postural need (as per prescription) and numerous pathologies throughout the whole life cycle of use.

In order to maintain its features it is necessary to periodically check the conservation and wear of all the parts composing the device as well as all hook and loops stripes. The positioning system seen the sizes and weight results easily transportable.

www.versainserto.com

1 MAN/WOMAN ADAPTABILITY

Finally a solution that makes the difference with a fully personalised fit.

2 ADAPTABILITY TO THE INDIVIDUAL MORFOLOGY AND ANATOMY

Unlimited possibilities to realize the customised postural solution.

3 PTS PELVIC TOTAL SUPPORT

A perfect fit to anatomy of the pelvis and related body segments, encouraging an active response of the musculo-skeletal system.

The distribution of the loads provided by the PTS reduces the mechanical stress on the deep tissues and significantly aids the prevention of decubitus ulcers.

4 FREE SHAPES PADDING

The innovative padding ensures maximum comfort and reliable protection.







MANUFACTURED BY:



Via Montagna Z.I. - Lotto 41
72023 Mesagne (BR) - ITALY
Tel. +(39) 0831 777840
Fax. +(39) 0831 730739
Email: sales@promedicare.it
www.versainserto.com
www.promedicare.eu

CERTIFIED COMPANY WITH QUALITY
MANAGEMENT SYSTEM ISO 9001 / EN ISO 13485
N° 279412 / N° 279411

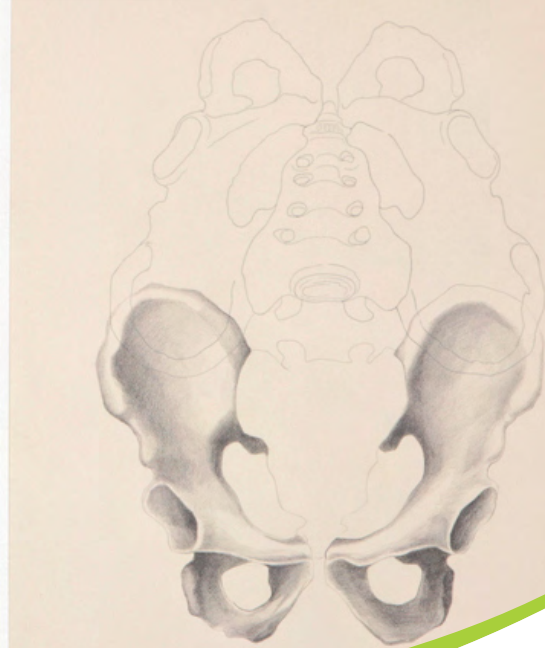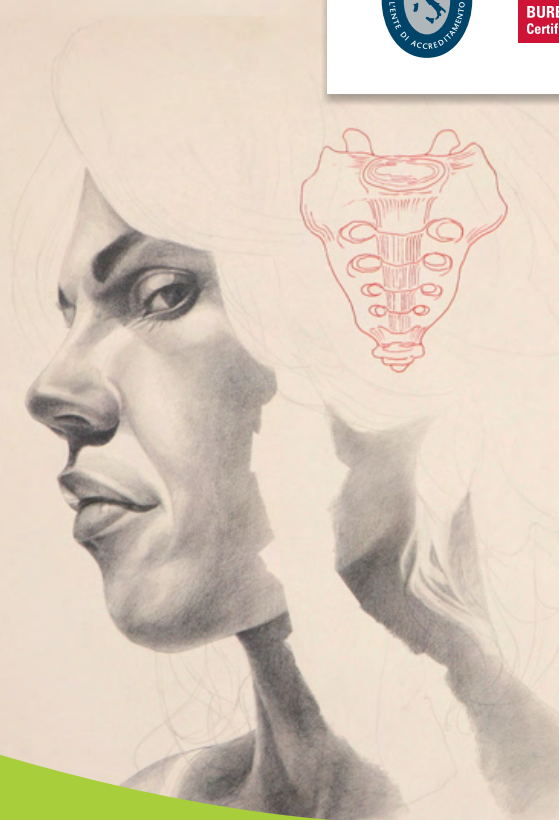

DISTRIBUTED BY:



CE 2 YEAR WARRANTY

www.versainserto.com