A NEW SOLUTION THE TRUNK TOTAL SUPPORT

A NEW SOLUTION TTS TRUNK TOTAL SUPPORT

A backrest without mounting hardware and shell





- FINALLY A BACKREST WITHOUT MOUNTING HARDWARE AND ALUMINIUM SHELL
- LIGHTWEIGHT AND COMFORTABLE
- FLEXIBLE TO MAXIMISE MOVEMENTS
- SUPPORT TO THE TRUNK IMMEDIATELY ADJUSTABLE WITHOUT TOOLS







CLINICAL INFORMATION

- Active users who need support in the lower part of the trunk to complete the stabilizing effects of the cushion.
- Mild instability of the lumbar and low thoracic parts of the spine.
- Users with an excellent control of the trunk.



TECHNICAL FEATURES

- Structural shaped shell in polyethylene closed cells.
- Three layers padding composed by a self-modelling polymeric foam, reticulated foam for ventilation and polyurethane open celled foam with low spring back action.
- Microclimatic ventilated cover.
- Multi-tension strap system for mounting on wheelchairs.

SPARE PARTS

Cover




For partially sighted users this document can be visualised in PDF format (with the possibility to enlarge the text) on the web site www.versainserto.com



FACILE BASIC



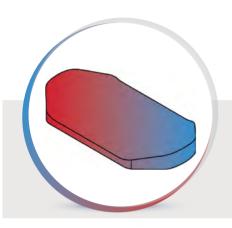

Inserto Facile Basic

IMMEDIATELY ADJUSTABLE WITHOUT TOOLS




MICROCLIMATIC VENTILATED COVER

° Each postural kit is ordered based upon the size of mobility/wheelchair device.

INSERTO FACILE BASIC SIZES TABLE										SIZES (cm.)		
MODEL	3415	3421	3618	3624	4018	4024	4418	4424	4820	4826		
WIDTH	34	34	36	36	40	40	44	44	48	48		
HEIGHT	15	21	18	24	18	24	18	24	20	26		
MOUNTING ON WHEELCHAIRS (OUTSIDE TUBES)	from 34 up to 36	from 34 up to 36	from 36 up to 38	from 36 up to 38	from 40 up to 42	from 40 up to 42	from 44 up to 46	from 44 up to 46	from 48 up to 50	from 48 up to 50		

EFFECTIVE DIMENSIONS TABLE											
SHELL EFFECTIVE WIDTH	28	28	30	30	34	34	38	38	42	42	
SHELL HEIGHT	15	21	18	24	18	24	18	24	20	26	

CONVERSION TABLE	INCHES	10	12	14	15	15,5	16	16,5	17,5	17,5	19	19,5
	CM.	25	30	36	38	39	40	42	44	45	48	50

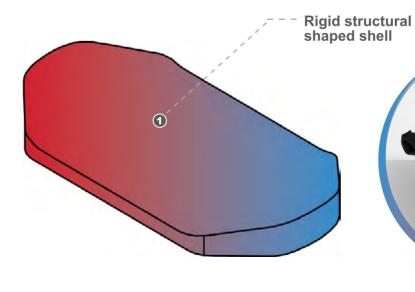





DISTRIBUTED BY:

CE 2 YEAR WARRANTY

www.versainserto.com