



Strong for **passion** 



















#### ADACTA LIGA GENERAL DESCRIPTION

Adacta Liga is the frame for positioning systems. It is made in alloy of anodised aluminum to guarantee lightness and durability.

The absence of welds makes it very safe.


The ADACTA LIGA frame, thanks to its modularity, presents a variety of configurations able to follow the pathological evolution of the user, the growth and consequential postural adjustments.

ADACTA LIGA can be set to the user's needs by acting on:


- The displacement of the upper frame (which houses the positioning system) with respect to the under frame
- The displacement of the rear wheels by guarantying the safety conditions.
   (as reported in the Instruction Manual)

#### ADACTA LIGA FRAMES FOR POSITIONING SYSTEMS

# "A CHOICE OF MODULARITY, VERSATILITY AND LIGHTNESS"

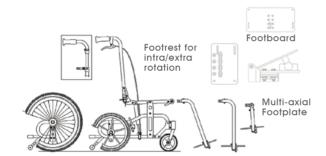

#### "A CHOICE OF SECURITY"

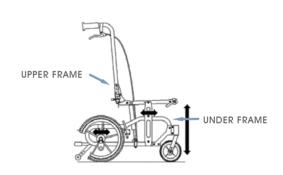



The most STABLE and RESISTANT between the frames for positioning systems



For transportation option it is strongly recommended the use of a headrest









### DIMENSIONS

SIZE (mm)	WIDTH	DEPTH
US	300	380
XXS	350	380
XS	350	420
XS1	380	460
S	410	460
M	410	500
M1	440	500
L	470	460
XL	470	500



Technical Features

- Supporting Frame Adacta Liga rigid or folding
- Continuous tilt in space from 0° up to 30° and on request -5° (by gas spring)
- Mechanical reclining (3°-9°-15°-22°)
- Folding down push handles
- Displacement of the upper frame with respect to the under frame
- Continuous displacement of the rear wheels plates
- Anti tip tubes
- Adjustable floor-to-seat-rails height (from 410 mm up to 490 mm.)

## OPTIONAL PARTS

Technical features

- Castors available:
- 125 mm. poly
- 150 mm. poly or pneumatic
- 175 mm. poly

- Quick release Rear Wheels available:
- 300 mm. poly or pneumatic
- 500 mm. poly or pneumatic with handrims
- 600 mm. pneumatic with handrims
- · Locks:
- With drum brakes
- With handle
- With drum brakes + handle
- Knee Angle:
- 90° or 75°
- Elevating with teethed mechanism
- Adduction/abduction adjustments
- Vertical and telescopic height adjustable
- Footplate:
- Separate footplates
- Separate footplates with intra-extra rotation slots
- Separate footplates with multi-axial mechanism
- Footboard
- Push Handles:
- Telescopic and rotating
- Transportation kit









### ACCESSORIES

- Rigid wooden base
- · Height adjustable flip back armrest
- Kit spondina detachable and multi-adjustable with multi-adjustable armrests
- Versa hip guides
- Versa arm supports
- Separate Calf support

- Calf panel
- Leg cradle
- Padded calf strap
- Positioning ankle support
- Versa Fixatis pelvic strap
- Versa Fixatis abduction harness
- Trays

The various range of accessories allows to set the system to user's needs.

### (HEIGHT (floor-to-seat-rails)

Different heights available (basing upon the wheels chosen) 410mm,430mm,450mm,470mm,490mm.

COLOR

ADACTA LIGA is available in 5 different colours







TELESCOPIC AND ROTATING PUSH HANDLES



SEAT/BACK MECHANICAL RECLINING WITH QUICK RELEASE PIN



LEVER FOR TILT MECHANISM



ANTI TIP TUBES



MULTI-ADJUSTABLE SPONDINA WITH ARM SUPPORT



RIGID WOODEN BASE



STABILISING CROSS BRACES



FOOTREST

## **TECHNICAL DIMENSIONS**

WIDT	H (mm)	US	XXS	XS	XS1	S	М	M1	L	XL
E :::	standard	480	530	530	560	590	590	620	650	650
300mm. WHEEL	w/drum brakes	490	540	540	570	600	600	630	660	660
یے غ	standard	500	550	550	580	610	610	640	670	670
500mm. WHEEL	w/drum brakes	530	580	580	610	640	640	670	700	700
600mm. WHEEL	standard	/	550	550	580	610	610	640	670	670

LOAD (kg)	FOLDING	RIGID
US	50	50
XXS	50	50
XS	50	50
XS1	75	75
S	75	75
M	75	100
M1	75	100
L	75	100
XL	75	100

FOLDING	RIGID		
12,2	11,6		
12,4	11,8		
12,8	12,2		
12,9	12,4		
13	12,5		
14	13,5		
14,5	14		
13,5	13		
15	14,5		
	12,2 12,4 12,8 12,9 13 14 14,5		



# **TECHNICAL DIMENSIONS**

C	OVERALL LENGTH (mm)		US	XXS	XS	XS1	S	M	M1	L	XL
300mm WHEEL	ELEVATING HANGERS	L	700	700	730	770	770	810	810	770	810
	75° HANGERS	L	680	680	710	730	730	760	760	730	760
		L.min	880	880	910	970	970	1000	1000	970	1000
등 뽀		L.max	990	990	1000	1040	1040	1070	1070	1040	1070
30 ×		L	680	680	710	730	730	760	760	730	760
	90° HANGERS	L.min	880	880	910	930	930	960	960	930	960
		L.max	940	960	990	1000	1000	1030	1030	1000	1030
	ELEVATING HANGERS	L	780	780	810	850	850	890	890	850	890
_	75° HANGERS	L	750	750	780	800	800	820	820	800	820
		L.min	950	950	980	1030	1030	1060	1060	1030	1060
500mm WHEEL		L.max	1030	1030	1070	1130	1130	1160	1160	1130	1160
50 W	90° HANGERS	L	750	750	780	800	800	820	820	800	820
		L.min	950	950	980	1000	1000	1030	1030	1000	1030
		L.max	970	970	1000	1060	1060	1090	1090	1060	1090
600mm WHEEL	ELEVATING HANGERS	L	/	830	840	880	880	920	920	880	920
	75° HANGERS	L	/	800	830	850	850	870	870	850	870
		L.min	/	1010	1030	1080	1080	1110	1110	1080	1110
		L.max	/	1090	1120	1160	1160	1200	1200	1160	1200
09 >	90° HANGERS	L	1	800	830	850	850	870	870	850	870
		L.min	/	1010	1030	1050	1050	1180	1180	1050	1180
		L.max	/	1050	1080	1120	1120	1160	1160	1120	1160

The length dimensions are subjected to variations due to the adjustment of the wheels plates. The dimensions Lmax and Lmin, referred to the elevating hangers, have not been reported because they depend on the inclination degree.

Graphic: www.huma.it



MANUFACTURED BY:



Via Montagna Z.I. Lotto 41 c.a.p. 72023 Mesagne (Br) - ITALY Tel. (+39) 0831.777840 Fax. (+39) 0831. 730739 Sito web: www.promedicare.it e mail: info@promedicare.it

### DISTRIBUTED BY:



The contents of this catalogue are property of Pro Medicare s.r.l. any improper or anauthorised use is strictly forbidden.